

Quality Resource Guide

Intraoral Bitewing Radiographic Technique

Author Acknowledgements

SALLY M. MAURIELLO, RDH BS EDD

Professor, Department of Dental Ecology

DAN SHUGARS, DDS PHD MPH

Professor Emeritus

University of North Carolina School of Dentistry
Chapel Hill, NC

Drs. Mauriello and Shugars have no relevant relationships to disclose.

The following commentary highlights fundamental and commonly accepted practices on the subject matter. The information is intended as a general overview and is for educational purposes only. This information does not constitute legal advice, which can only be provided by an attorney.

© Metropolitan Life Insurance Company, New York, NY.
All materials subject to this copyright may be photocopied for the noncommercial purpose of scientific or educational advancement.

Originally published July 2002. Updated and revised October 2005, May 2008, May 2011, December 2014 and December 2017. Expiration date: December 2020.

The content of this Guide is subject to change as new scientific information becomes available.

Educational Objectives

Following this unit of instruction, the dental team member should be able to:

1. Describe the characteristics of a diagnostic bitewing radiograph.
2. Recognize the guidelines for proper receptor placement and beam orientation.
3. Recognize the common errors found in bitewing projections and their cause.
4. Describe the ways in which common errors can be avoided.

Introduction

Bitewing radiographs are the most common intraoral radiographic tool used for diagnosis in dentistry. It is estimated that bitewings account for nearly 60% of all dental radiographic surveys taken. When combined with thorough clinical examinations, bitewings are a valuable diagnostic tool. They are excellent aids for the identification of certain types of pathology, defective restorations, dental caries, and periodontal disease. However, like most radiographic procedures, bitewings are technique-sensitive. In an attempt to prevent or reduce the frequency of technique-related problems, this module will describe the characteristics of good diagnostic bitewings, present three common technique errors, and describe methods for solving technique errors.

MetLife designates this activity for **1.0 continuing education credits** for the review of this Quality Resource Guide and successful completion of the post test.

ADA CERP® | Continuing Education Recognition Program

MetLife is an ADA CERP Recognized Provider.

Accepted Program Provider FAGD/MAGD Credit **11/01/16 - 12/31/20**.

ADA CERP is a service of the American Dental Association to assist dental professionals in identifying quality providers of continuing dental education.

ADA CERP does not approve or endorse individual courses or instructors, nor does it imply acceptance of credit hours by boards of dentistry.

Address comments or questions to:

DentalQuality@metlife.com

MetLife Dental Continuing Education
501 US Hwy 22, Area 3D-309B
Bridgewater, NJ 08807

Cancellation/Refund Policy:

Any participant who is not 100% satisfied with this course can request a full refund by contacting us.

Concerns or complaints about a CE provider may be directed to the provider or to ADA CERP at www.ada.org/goto/cerp.

Bitewing Devices

There are primarily two common devices used to expose bitewing projections: bitewing tabs and XCP® bitewing holders. The technique and strengths/weaknesses are discussed for both.

Bitewing Tabs


The bitewing tab is attached to the active side of the receptor cover and the circular collimator aligned to completely irradiate the receptor. The collimator is preset at a positive 10 degrees angulation. The disadvantages of this technique are: the inability to use rectangular collimation and the freehand alignment of the beam to the receptor.

Beam-Alignment Device (XCP®)

Use of a beam-alignment receptor holder allows for the use of a rectangular collimator which decreases the radiation dose to the patient. The receptor is placed in the holder with the active side facing the collimator. Once the receptor is placed orally, as described in this module, the rectangular collimator is aligned to the positioning collimator. The positioning collimator automatically aligns the beam at 90 degrees to the receptor. The one limitation to using the beam alignment device is the potential to sacrifice horizontal alignment resulting in interproximal overlap. To address this limitation, place the holder and then rotate the aiming ring horizontally to open the interproximal contact.

Characteristics of Diagnostic Bitewing Projections

The radiographic projection best suited for the interpretation of dental caries and periodontal disease in the posterior area of the mouth is the bitewing radiograph. This is primarily due to the vertical angulation used in exposing the projection. In contrast, in a periapical projection

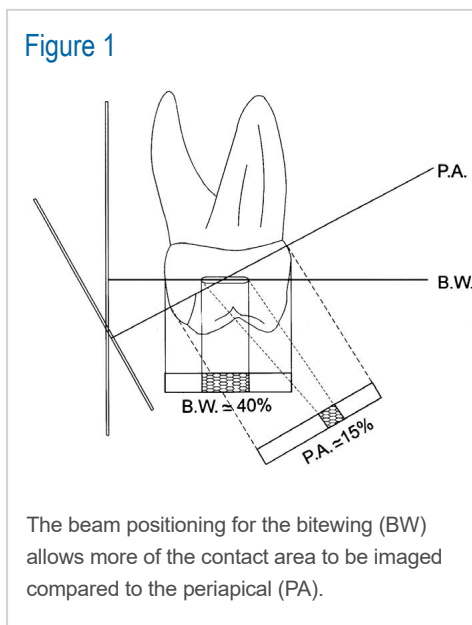

the vertical angulation can distort the caries lesion or superimpose other structures to mask or hide a lesion or distort bone topography.

Figure 1 depicts the projection geometry principles for bitewing and periapical radiographs. With the angulation of a bitewing projection, approximately 40% of the proximal surface is imaged compared to only 15% in a periapical projection.

Because radiographs can detect proximal caries lesions and alveolar bone not easily seen in a clinical examination, it is essential that the radiograph clearly display interproximal surfaces of the teeth and the crestal ridge of the adjacent alveolar bone. For a caries lesion to be detected radiographically, the mineral content of the tissue must change or demineralize 30-50%. As a result, a lesion will often tend to be larger clinically than that displayed by the radiographic image. Thus, an image with adjacent interproximal surfaces not overlapped will enhance the ability to detect lesions early, which in turn impacts management considerations.

Problem-Solving Techniques

Receptor Placement

Packet placement is one of the most critical elements in exposing a diagnostic radiograph. It is also one of the most common sources of error. Because of the contour of the dental arch, it is often difficult to place the packet and angle the beam to create an open proximal contact between those teeth. There are two easy steps that ensure proper placement of the receptor. First, place the receptor as far forward as possible. Second, angle the receptor so it is positioned behind the mandibular lateral on the opposite side of the mouth. Remember, for improved patient comfort and ease in positioning the packet, place the packet towards the tongue or midline of the floor of the mouth by pushing towards the tongue. This also allows for proper positioning behind the mandibular lateral on the opposite side. This technique is also effective when tori are present.

Beam Orientation

Proper orientation of the central beam is essential to obtaining diagnostic projections. Misalignment of the x-ray beam through the interproximal spaces can cause overlap of the proximal areas. When this occurs, the SLOB rule (Same Lingual Opposite Buccal) can be used to determine the angulation that caused the error and thus correct the error on the retake projection. To identify the cause of the horizontal overlap, the clinician can use the following rule. First, identify the lingual cusp on either maxillary premolar on the receptor. This would be the shorter cusp. Then, determine whether the lingual cusp is mesial or distal to the facial cusp (the longer cusp). Using the SLOB Rule when viewing a radiograph, the tubehead would be positioned too mesial if the lingual cusp is mesial to the facial cusp and too distal if the lingual cusp is distal to the facial cusp. Thus, the horizontal overlap error can be corrected by moving or repositioning the tubehead in the opposite direction.

Notice that the radiographs shown in **Figures 2** and **3** meet the four criteria of diagnostic premolar and molar bitewings. These criteria are:

- For the premolar projection, the distal of the canine (**Figure 2**) and in the molar projection no more than the distal of the mandibular second premolar should be imaged (**Figure 3**).
- The interproximal spaces should be open with no horizontal overlap between the maxillary premolars on the premolar projection. For the molar projection, there should be no horizontal overlap between the maxillary 1st and 2nd molars.


- The alveolar crestal bone should be imaged with about the same amount of bone in both arches displayed on the receptor.

- For the molar projection, all of the erupted terminal molar should be displayed with the retromolar area imaged.

Figure 2



Figure 3

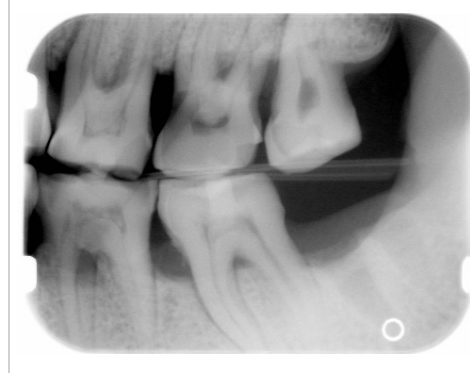

Common Errors

Figure 4



Problem:

This projection did not image the mesial of the first premolars onto the receptor due to improper packet placement in the mouth.

Solution:

This error can be corrected by placing the anterior edge of the receptor behind the mandibular lateral on the opposite side of the arch.

Figure 5



Problem:

The mesial of the premolars is missing on this projection due to improper horizontal angulation. The tubehead was too distal, thus projecting the mesial of the premolars off of the receptor.

Solution:

This error can be corrected by moving the tubehead towards the mesial.

Figure 6



Problem:

The third molar is not imaged on the receptor due to improper receptor placement.

Solution:

Align the front edge of the receptor so that it is no further forward than the mesial of the mandibular first molar. If the tooth is still not imaged on the receptor, then a third molar disto-oblique projection should be exposed.

Common Errors (continued)

Figure 7



Problem:

The interproximal areas are overlapped, thus obscuring any carious lesions. In this example, the tubehead is directed too mesial, determined by the lingual cusp (shorter cusp) positioned mesial to the buccal cusp (longer cusp).

Solution:

Rotate the tubehead toward the distal.

Figure 8



Problem:

The level of alveolar bone imaged in the projection is unequal in the maxillary and mandibular arches. This error could be due to the receptor-holding device not centered on the receptor, or incorrect vertical angulation positioned too steep (tubehead positioned at an angle greater than +10).

Solution:

Reposition the receptor so that it is centered in the holder. The vertical angulation should be positioned at +10 degrees (beam is angled down). If the patient exhibits periodontal pocketing of 5 mm or more, then a vertical bitewing should be exposed to ensure that the crestal bone is imaged on the projection.

Summary

Bitewing projections play an important role in aiding the identification and treatment of dental disease. Therefore, it is imperative that good diagnostic bitewing radiographs be exposed. This guide describes the desired characteristics of posterior bitewings, problem-solving techniques and common problems with associated solutions. It was written based upon existing guidelines, however, organizations such as NCRP often provide periodic updates. Readers are encouraged to source the appropriate agency or organization for the most current recommendations.

Suggested Reading


Oral Radiology: Principles and Interpretation, Chap. 7, 14 and 16 by White and Pharoah, 7th Ed. Mosby, 2014

Dental Radiographic Examinations: Recommendations for Patient Selection and Limiting Radiation Exposure. American Dental Association, Council on Scientific Affairs and US Department of Health and Human Services, Public Health Service, Food and Drug Administration, 2012.

POST-TEST

Internet Users: This page is intended to assist you in fast and accurate testing when completing the "Online Exam." We suggest reviewing the questions and then circling your answers on this page prior to completing the online exam.

(1.0 CE Credit Contact Hour) Please circle the correct answer. 70% equals passing grade.

1. **The use of beam alignment devices for bitewing radiography is recommended so that:**
 - a. rectangular collimation can be used to reduce patient dose
 - b. the technology improves horizontal angulation
 - c. cone cuts don't occur
 - d. the technique is more comfortable to the patient
2. **A molar bitewing should image no more than the distal of the mandibular second premolar.**
 - a. True
 - b. False
3. **The correct vertical angulation to use with a bitewing is:**
 - a. +10 degrees
 - b. +20 degrees
 - c. -10 degrees
 - d. -20 degrees
4. **Vertical bitewings are indicated when there is ___ mm of periodontal pocketing.**
 - a. 3
 - b. 4
 - c. 5
 - d. All of the Above
5. **The use of bitewing tabs requires all EXCEPT:**
 - a. freehand alignment of the beam to the receptor
 - b. the use of rectangular collimation
 - c. the use of circular collimation
 - d. presetting the vertical angulation to +10 degrees
6. **A bitewing radiograph that exhibits unequal amounts of alveolar bone in each of the arches is most likely a result of:**
 - a. incorrect vertical angulation
 - b. incorrect horizontal angulation
 - c. Both a and b
 - d. None of the Above
7. **A molar bitewing radiograph that displays a partial third molar, most likely is a result of:**
 - a. incorrect vertical angulation
 - b. incorrect horizontal angulation
 - c. incorrect packet placement
 - d. None of the Above
8. **Bitewings are used to aid in the identification of:**
 - a. defective restorations
 - b. caries lesions
 - c. periodontal disease
 - d. All of the Above
 - e. None of the Above
9. **If the mesial of the first premolar is not imaged on the projection and there is no horizontal overlap, then the most likely cause for the missing mesial of the premolar is:**
 - a. incorrect vertical angulation
 - b. incorrect horizontal angulation
 - c. incorrect packet placement
 - d. None of the Above
10. **Using the SLOB rule, determine the cause of the horizontal overlap in the maxillary premolar area in this radiograph.**

 - a. the tubehead is too mesial
 - b. the tubehead is too distal
 - c. the tubehead is positioned correctly, the receptor is turned
 - d. None of the Above

Registration/Certification Information (Necessary for proper certification)

Name (Last, First, Middle Initial): _____

Street Address: _____ PLEASE PRINT CLEARLY Suite/Apt. Number _____

City: _____ State: _____ Zip: _____

Telephone: _____ Fax: _____

Date of Birth: _____ Email: _____

State(s) of Licensure: _____ License Number(s): _____

Preferred Dentist Program ID Number: _____ Check Box If Not A PDP Member

AGD Mastership: Yes No

AGD Fellowship: Yes No Date: _____

Please Check One: General Practitioner Specialist Dental Hygienist Other

**FOR
OFFICE
USE
ONLY**

Evaluation - Intraoral Bitewing Radiographic Technique 6th Edition

Providing dentists with the opportunity for continuing dental education is an essential part of MetLife's commitment to helping dentists improve the oral health of their patients through education. You can help in this effort by providing feedback regarding the continuing education offering you have just completed.

Please respond to the statements below by checking the appropriate box, using the scale on the right.

1 = POOR

5 = Excellent

	1	2	3	4	5	
1. How well did this course meet its stated educational objectives?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
2. How would you rate the quality of the content?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
3. Please rate the effectiveness of the author.	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
4. Please rate the written materials and visual aids used.	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
5. The use of evidence-based dentistry on the topic when applicable.	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/> N/A
6. How relevant was the course material to your practice?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
7. The extent to which the course enhanced your current knowledge or skill?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
8. The level to which your personal objectives were satisfied.	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
9. Please rate the administrative arrangements for this course.	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	

10. How likely are you to recommend MetLife's CE program to a friend or colleague? *(please circle one number below:)*

10
9
8
7
6
5
4
3
2
1
0

extremely likely
neutral
not likely at all

What is the primary reason for your 0-10 recommendation rating above?

11. Please identify future topics that you would like to see:

Thank you for your time and feedback.



To complete the program traditionally, please mail your post test and registration/evaluation form to:
MetLife Dental Quality Initiatives Program | 501 US Highway 22 | Bridgewater, NJ 08807